



BIOPRESS™ CULTURE PLATES

Flexible bottomed culture plate for use with the Flexcell Compression system.

- Each well contains a stationary platen and piston between which 3-D tissue samples or cell seeded constructs are deformed using the FX-5000 Compression system.
- Ring foam holders keep samples confined to central region of each well.
- All plates come pre-sterilized.
- Available individually or in cases of 40 plates.




Ordering Information	
<i>Catalog Number</i>	<i>Product/Item</i>
BF-3000C	BioPress Compression Culture Plate

STAGEPRESSER™ MEMBRANES

Flexible growth surface for use with the StagePresser device.

- Compatible with the Flexcell Compression system.
- Each membrane contains an acrylic piston which is compressed against a stationary platen on the StagePresser device to deform 3-D tissue samples or cell-seeded constructs.
- Ring foam holders keep samples confined to central region of each membrane.
- Each membrane comes individually packaged in a sterile cell culture dish.

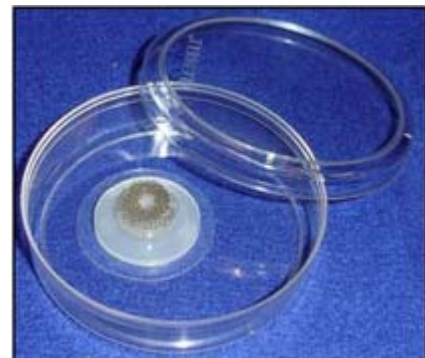

Ordering Information	
<i>Catalog Number</i>	<i>Product/Item</i>
SPM-3000	StagePresser Membrane (Piston and Foam Retainer)