

## Laboratory equipment suitable for gel and blot applications

### Shaker "Belly Dancer"



- Suitable for gel and blot applications
- Continuous gentle or thorough shaking
- Unique three dimensional motion guarantees accurate mixing action and offers the possibility to stack blotting trays
- Durable construction, low noise level
- Ambient temperature operation: 0 °C to 40 °C
- Can be used in incubators and cool rooms
- Speed control: up to 100 rpm
- Maximum load: 9.1 kg
- Platform size: 300 x 300 mm
- Dimensions (w x d x h): 381 x 381 x 230 mm
- Additional platform can be added (see picture)

Cat.-No.	Description
BDRLS0006	"Belly Dancer", with finger mat for tubes
BDRLS0004	"Belly Dancer", with additional platform (approximately 300 x 300 mm)
BDRLS0005	"Belly Dancer", with extended platform (approximately 450 x 450 mm)

### Shaker "Belly Button"



- Suitable for gel and blot applications
- Very compact version of the "Belly Dancer" (see above), fits also into smaller incubators
- Platform tilt is continuously adjustable from level to 8 °
- Ambient temperature operation: 0 °C to 60 °C
- Can be used in incubators and cool rooms
- Speed control: 6 - 100 rpm
- Maximum load: 2.3 kg
- Dimensions (w x d x h): 203 x 203 x 121 mm
- Optionally available: non-skid mat for shaking tubes (see picture)

Cat.-No.	Description
BBUAAUV1S	"Belly Button"
ACCBB0001	Non-skid mat for tubes

### Roller Mixers



- Space saving low profile design
- Available with 6 or 11 rolls (6 rolls model can be extended to 11 rolls)
- Suitable for culture and hybridization bottles or tubes due to removable rolls, dimension of one roll (Ø x L): 29 x 254 mm
- Can be used with 3 big culture bottles (850 cm<sup>2</sup>), 4 hybridization bottles or with up to 30 tubes (10 x 75 mm), if used with 11 rolls
- Speed control: 2 - 38 rpm
- Maximum load: 4 liter liquid
- Ambient temperature operation: 0 °C to 60 °C
- Dimensions (w x d x h): 362 x 381 x 60.5 mm

Cat.-No.	Description
ROLAAUV1S	Roller mixer with 6 rolls
ROLLS0001	Roller mixer with 11 rolls
ACCRLSR10	Single roll

## Shaking Water Bath



- Combination of modified "Belly Dancer" (see front page) and hybridization water bath
- Units can also be used separately
- Separation of shaker and water bath offers the possibility for easy cleaning and decontamination when working with radioactive material
- PETG gable bath lid (thickness: 5 mm) against evaporation and beta radiation
- PID microprocessor heat control, safety thermostat
- Temperature range: 30 °C to 75 °C
- Quick heating, 3 °C rise per minute
- Excellent temperature control and stability:  $\pm 0.2$  °C

Cat.-No.      Description  
 BWBAA220S    **Shaking Water Bath**

## Controller "Blot Washer"




**Belly Dancer**  
 (has to be ordered separately)

- The "Blot Washer" is a controller with integrated pump, ideal for the automation of multiple washing steps for Western Blots
- Up to 10 washing programmes can be created and saved
- Programmable washing time from seconds to 99 hours
- Simultaneous processing of up to four containers with different sizes
- Liquid dispensing rate: approx. 2 ml/second
- Liquid aspirating rate: approx. 5 ml/second
- Delivered with 2 bottles of 4 liter and 2 blot trays
- Can be combined with the "Belly Dancer" (see front page, not included, has to be ordered separately)
- Accessories for "Blot-Washer" ( tubing, clamps, clips, trays, etc.) on request

Cat.-No.      Description  
 BLWAA220S    **"Blot Washer"**, controller with integrated pump

## "Hi/Lo Profile Rocker"

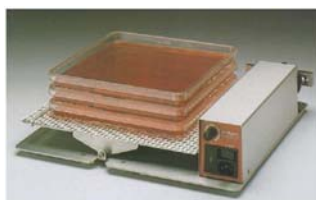

- Space saving low profile design, stainless steel
- Autoclavable (without motor)
- Gentle to vigorous shaking motion
- Speed control: 0 - 22 motions/minute
- Platform tilt is continuously adjustable from: 2 ° - 11 °
- Ambient temperature operation: 0 °C to 60 °C
- Can be used in incubators and cool rooms
- Maximum load: 9.1 kg
- Platform with non-skid mat: 305 x 305 mm
- Option: extended capacity with up to 3 additional platforms
- Distance between platforms (if more than one is used): 114 mm
- Dimensions (w x d x h): 349 x 400 x 44.5 mm (without additional platforms)

Cat.-No.      Description  
 ROCAA220S    **"Hi/Lo Profile Rocker"**  
 ROCLS0004    **"Hi/Lo Profile Rocker"**, with 1 additional platform  
 ROCLS0005    **"Hi/Lo Profile Rocker"**, with 2 additional platforms  
 ROCLS0006    **"Hi/Lo Profile Rocker"**, with 3 additional platforms

**Please note that in Europe all equipment listed above is not delivered with a flat cable!**