

COMPACT POWER. CUSTOM FLEXIBILITY.

Versatile speed, swappable trays, and timed control, all in one durable, space-saving shaker.

One size fits most.
CUSTOM fits the rest.

OPTIONS
Standard
Mini


TRAY OPTIONS

10" x 10" (25 x 25 cm) Tray
20" X 20" (51 x 51 cm) Tray
Mini-Shaker Tray 8-Position
Mini-Shaker Tray 32-Position

AVAILABLE VOLTAGE

100V
115V
230V

MAXIMUM LOAD
24 lbs (11 kg)

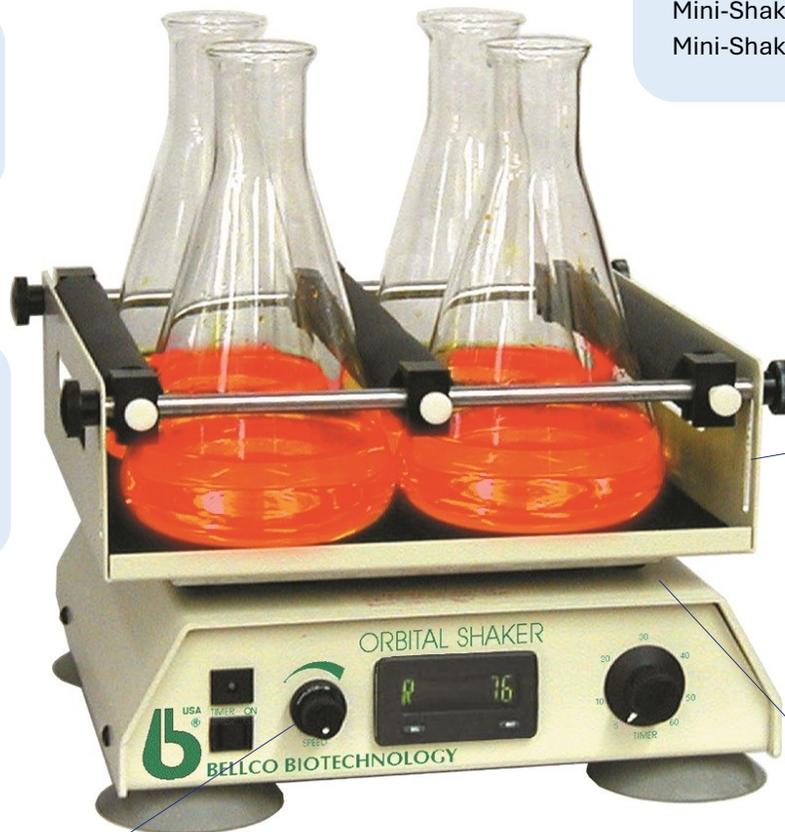

SPEED RANGE

0 – 350 rpm
40 – 1550 rpm

ORBIT

1" (25 mm)
1/16" (1.6 mm) - Mini

MADE IN 
USA 

*Most Bellco products can be configured to suit your unique needs.

DIGITAL ORBITAL SHAKER

The Bellco Digital Shaker is a compact, heavy-duty unit for agitating gels and assay media, perfect for staining and destaining. It offers variable speeds from 30–350 rpm and supports a range of vessels with interchangeable trays, including flask clamp and peg-and-grommet styles. A timer allows automatic shutoff from 2–120 minutes.



Specifications:

- Speed Range: 30 – 350 rpm +/- 2%
- Circular Orbit: 1" (25 mm)
- Base Overall Dimensions:
12" x 14" x 4 1/2"
(13 x 28 x 28 cm)

DIGITAL MINI ORBITAL SHAKER

The Bellco Digital Mini-Orbital Shaker is a high speed shaker for general laboratory cell culture experimentation. The design incorporates an orbital platform that accepts a variety of optional trays. Peg and grommet design offers a simple method to install or change 4 versions of trays. It has a brushless DC motor designed for minimal maintenance.



Specifications:

- Speed Range: 40 – 1550 rpm
- Circular Orbit: 1/16" (1.6 mm)
- Base Overall Dimensions:
13 3/4" x 14 1/2" x 5 1/2"
(35 x 37 x 14 cm)



CONTACT US:

Tel.: +49 (0) 26 83 / 4 30 94
E-mail: info@dunnlab.de
Internet: www.dunnlab.de

